

ChM produces and distributes advanced medical solutions in 3 main divisions:



Comprehensive portfolio of products for traumatology, including systems designed for fractures fixation and deformities correction of extremities and pelvis.



Wide range of advanced solutions for cervical and thoraco-lumbar stabilization of spine, including pedicle screw systems for open and MIS procedures, various interbody devices and fixation plates.



Instruments and implants for cranio-maxillofacial surgeries, dedicated for fracture fixations, reconstructions, distractions and orthognathic surgeries.



NARROW STRAIGHT PLATES

3.7093.6xx



To bring medical solutions

Document No B/P-18
Review date P-002:19.08.2020

www.chm.eu

NARROW STRAIGHT PLATES



NARROW STRAIGHT PLATES

Design aligned with anatomy

- optimized anatomical profile
- facilitated anatomical reduction

Bottom undercuts of the shaft part

- limited bone-to-plate contact
- better blood circulation of periimplant tissues

Beveled tip

- easy percutaneous insertion

Compression hole




- oblong hole for plate positioning
- compression in two directions possible

K-wire holes

- provisional plate positioning

Multiple screw options

- non-locking and locking screw gives multiple configuration for individual cases

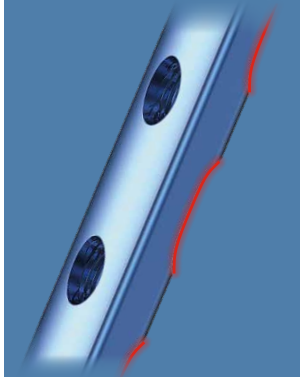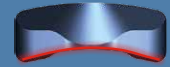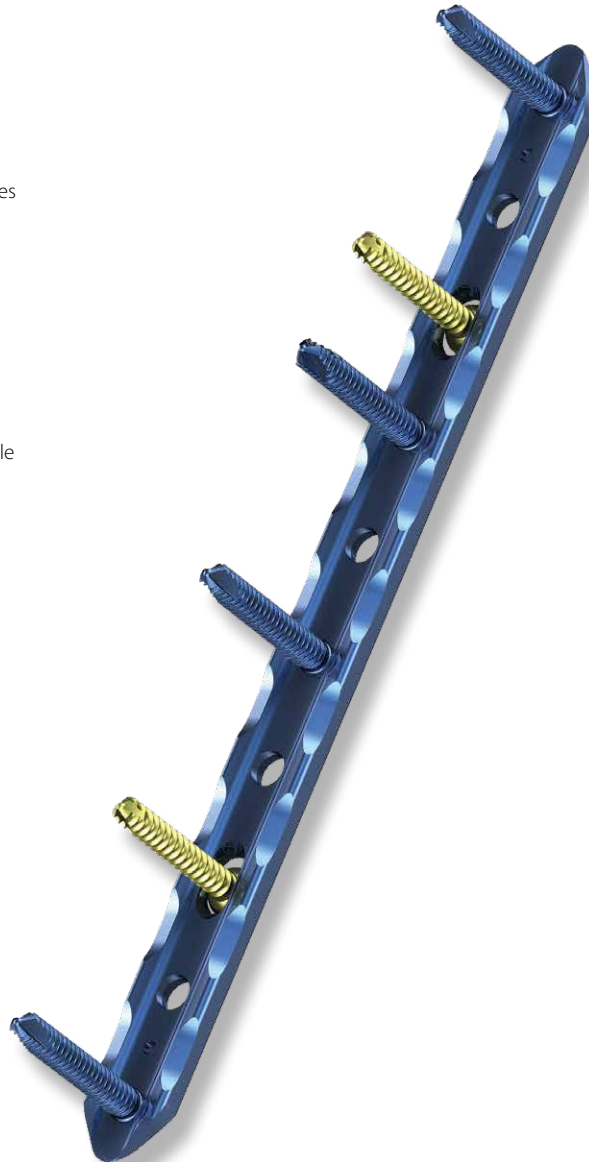

Design aligned with anatomy

Bottom undercuts of the shaft part

Compression/K-wire holes

Beveled tip